




# ZAPI LUXE™



## PRODUCT FEATURES

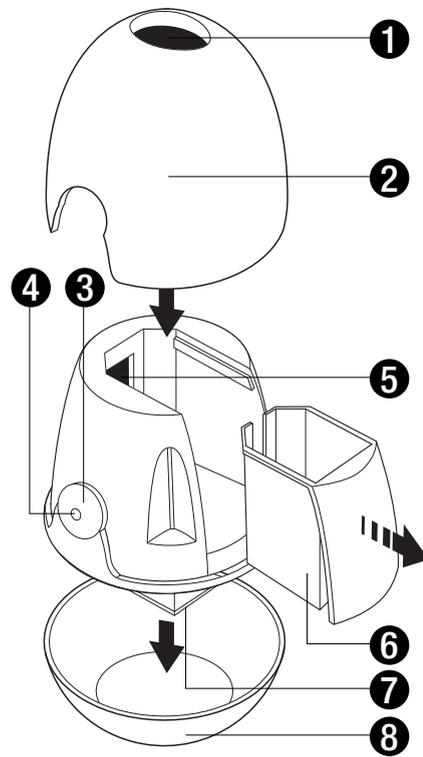

1. Toothbrush Insertion Slot
2. Removable Lid
3. Power Button
4. Low Battery Light
5. UV Sanitizing Light
6. Removable Drip Cup/  
Toothbrush Holder
7. Battery Compartment
8. Removable Weighted Base

## BOX CONTENTS

- 1 Zapi Luxe™
- 1 Toothbrush

## GETTING STARTED (Inserting Batteries)

1. With the unit turned upside down, rotate the bottom base counterclockwise to unlock then remove base. Open the battery cover and place 3 “AA” batteries in the slots as shown (**Diagram 1**).

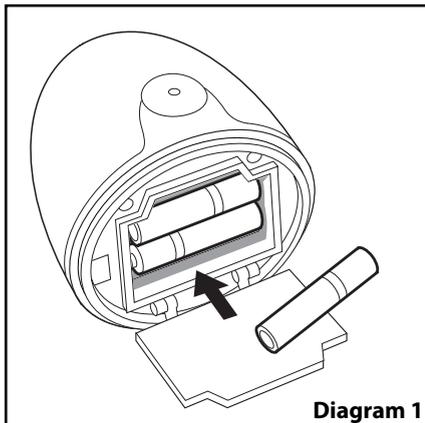


Diagram 1

2. Shut lid of battery compartment.

3. Align the bottom base to the top of the Zapi Luxe and turn clockwise until it snaps closed (**Diagram 2**).

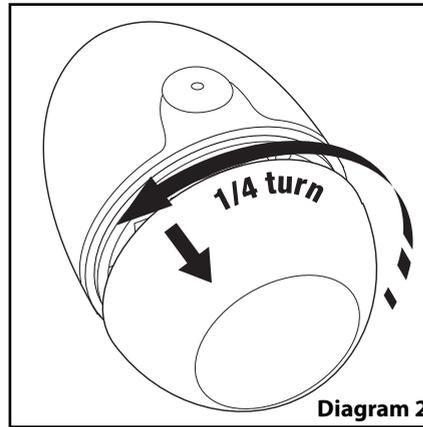


Diagram 2

## USING ZAPI LUXE™

1. With the power button of the Zapi Luxe facing you, insert your toothbrush or electric brush head with the bristles facing left into the top hole of your Zapi Luxe.
2. When the head of the toothbrush has been inserted into the Zapi Luxe gently click the handle into the rubber holder (**Diagram 3**).

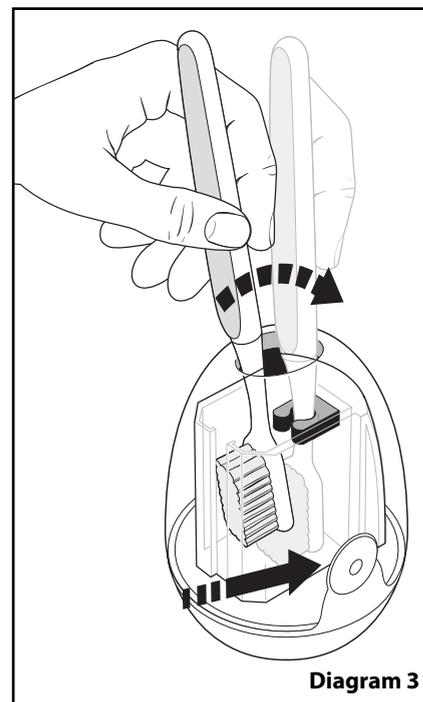


Diagram 3

3. Press the button on the front of the Zapi Luxe to activate the UV light. The blue light will blink slowly for approximately 6 minutes during the sanitizing cycle.

4. Once the cycle is complete, the toothbrush can remain in the Zapi Luxe until your next brushing.

5. Zapi Luxe has a low battery indicator light. When the light in the center of the power button begins to blink quickly it is time to change the batteries.

## Maintenance & Cleaning

1. Drip cup should be cleaned weekly. Depress sides of Zapi Luxe cover and pull up (**Diagram 4**). If you are having difficulty, insert a coin at the seam and pry up to loosen. Pull up the drip cup and clean with soap and warm water. Let air dry before reinserting into the Zapi Luxe.

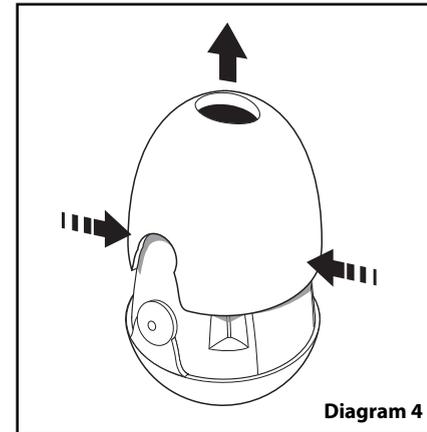


Diagram 4

**NOTE: Tap off any excess water from your toothbrush bristles before inserting it into the Zapi Luxe.**

2. The outside of the Zapi Luxe can be cleaned with a damp cloth.

## Troubleshooting

If the unit is not turning on, check to see if the batteries were inserted correctly.

If the batteries are inserted correctly and the unit still does not turn on, make sure that the lid is snapped closed.

Do not mix old and new batteries. Replace both when the unit is no longer functioning.

## UV Toothbrush Sanitizer

VZL100 Product Series